



User's Manual

Fenspiride ELISA Kit



DEIA-XYL29



96T



This product is for research use only and is not intended for diagnostic use.

For illustrative purposes only. To perform the assay the instructions for use provided with the kit have to be used.

Creative Diagnostics

 Address: 45-1 Ramsey Road, Shirley, NY 11967, USA

 Tel: 1-631-624-4882 (USA) 44-161-818-6441 (Europe)  Fax: 1-631-938-8221

 Email: info@creative-diagnostics.com  Web: www.creative-diagnostics.com

PRODUCT INFORMATION

Intended Use

For the determination of trace quantities of Fenspiride and/or other metabolites in human urine, blood, oral fluid.

General Description

CD's Fenspiride ELISA (Enzyme-Linked ImmunoSorbent Assay) test kit is a qualitative one-step kit designed for use as a screening device for the detection of drugs and/or their metabolites. The kit was designed for screening purposes and is intended for forensic use only. It is recommended that all suspect samples be confirmed by a quantitative method such as gas chromatography/mass spectrometry (GC/MS).

Principles of Testing

CD's test kit operates on the basis of competition between the drug or its metabolite in the sample and the drug-enzyme conjugate for a limited number of antibody binding sites. First, the sample or control is added to the microplate. Next, the diluted drug-enzyme conjugate is added and the mixture is incubated at room temperature. During this incubation, the drug in the sample or the drug-enzyme conjugate binds to antibody immobilized in the microplate wells. After incubation, the plate is washed 3 times to remove any unbound sample or drug-enzyme conjugate. The presence of bound drug-enzyme conjugate is recognized by the addition of K-Blue[®] Substrate (TMB). After a 30 minute substrate incubation, the reaction is halted with the addition of Red Stop Solution. The test can be read visually or with a microplate reader equipped with a 650 nm filter. The extent of color development is inversely proportional to the amount of drug in the sample or control. In other words, the absence of the drug in the sample will result in a dark blue/purple color, whereas the presence of the drug will result in light pink to no color development.

Reagents And Materials Provided

1. EIA Buffer: 40 mL (ready-to-use). Phosphate buffered saline solution with bovine serum and a preservative. Provided for dilution of drug-enzyme conjugate concentrate and samples.
2. Wash Buffer Concentrate (10X): 20 mL. Phosphate buffered saline solution with a surfactant. Dilute 10 fold with deionized or ultrapure water before use. Diluted wash buffer is used to wash all unbound conjugate and samples from the plate after the conjugate incubation.
3. K-Blue Substrate: 20 mL (ready-to-use). Stabilized 3, 3', 5, 5' Tetramethylbenzidine (TMB) plus Hydrogen Peroxide (H₂O₂) in a single bottle. It is used to develop the color in the wells after washing. Light sensitive.
4. Drug-Enzyme Conjugate: 2 x 100 µL. Lyophilized vials of Drug-horseradish peroxidase concentrate. Dilute 180X before use.
5. Antibody Coated Plate: A 96 well Costar plate, in strips of 8 break-away wells, coated with anti-drug antiserum. The plate is ready for use as is. Do not wash.
6. Red Stop Solution: 20 mL (ready-to-use). Non-acidic reagent used to stop the enzyme reaction.
7. Qualitative QC Positive Control: 750 µL provided (synthetic human urine). Do not dilute.

8. Qualitative QC Negative Control: 750 µL provided (synthetic human urine). Do not dilute.

Materials Required But Not Supplied

1. Deionized water.
2. Precision pipettes that range from 10 µL - 1000 µL and disposable tips.
3. Graduated cylinder to dilute and mix wash buffer.
4. Plate cover or plastic film to cover plate during incubation.
5. Clean glassware (i.e. test tubes) to dilute drug-enzyme conjugate.
6. Microplate reader with 650 nm filter if Red Stop is used, or a 450 nm filter if 1N HCl is used to stop the reaction. Note: It is not necessary to stop the reaction if reading immediately. Unstopped reactions should be read with a 650 nm filter.
7. Cut-off calibrator.

OPTIONAL MATERIALS

1. 1N HCl if Red Stop Solution is not used.
2. Microplate shaker.

Storage

This kit can be used until the expiration date on the label when stored refrigerated at 2-8°C. Store controls frozen if not used within 10 days. Note: Some kits require controls to be stored frozen immediately upon receipt. Reference kit label for details.

Specimen Collection And Preparation

Recommended minimum sample dilutions are listed below. These dilutions may change based on your laboratory's determination. All sample dilutions should be made in CD's EIA Buffer.

- a. Urine: A dilution with EIA Buffer may be necessary to reduce natural background as well as bring desired cutoff concentration within the assay range. Please contact your CD Representative for assistance.
- b. Whole blood: A dilution of 1:5 (i.e. 1 part sample to 4 parts provided EIA Buffer) is recommended. Please contact your CD Representative for assistance.
- c. Other Forensic sample types: Please contact your CD Representative for assistance.

Assay Procedure

PROCEDURAL NOTES

1. Desiccant bag must remain in foil pouch with unused strips. Keep zip lock pouch sealed when not in use to maintain a dry environment.
2. Use clean pipette tips for the buffer, drug-enzyme conjugate, controls and samples.
3. Before pipetting a reagent, rinse the pipette tip three times with that reagent.
4. When pipetting into the wells, DO NOT allow the pipette tip to touch the inside of the well or any of the



reagent already inside the well. This may result in cross contamination.

5. Controls and samples should be assayed in duplicate.
6. The drug-enzyme conjugate is most stable in its concentrated form. Dilute only the volume necessary for the amount of strips currently being used.
7. Before opening the drug-enzyme conjugate vial, tap the vial in an upright position to remove any liquid in the cap.
8. Before substrate addition, wipe the outside bottom of the wells with a lint-free wiper to remove dust and fingerprints.
9. Gently mix specimens and reagents before use. Avoid vigorous agitation.

TEST PROCEDURES

The following test procedures can be run manually or on an automated instrument. Please contact your CD representative for assistance with protocols for automated instruments.

1. To reconstitute the lyophilized conjugate, pipette 100 µL of deionized water into the vial. Allow the conjugate to rehydrate by slowly rolling the vial. Invert the vial (tapping on the cap) to reconstitute any dry material trapped in the cap. Bring vial to an upright position and again tap the cap to release the volume from the cap. Do not vortex or shake. Avoid excess foaming. Continue to invert until all solid material is reconstituted. Allow the conjugate concentrate to incubate for 15 minutes. The shelf life of the reconstituted conjugate is 7 days. The conjugate will reach full activity upon completion of this incubation. Label the vial with the date of reconstitution for future reference. Once vial is reconstituted, store at 2-8°C. Do not freeze. Note: If more concentrated conjugate is needed then is contained in the first conjugate vial, reconstitute and use the second vial. Do not use the separate contents of both vials in the same assay as some assay variability may result. If the contents of both vials are required of an assay, pool the concentrated conjugates. Use the expiration date of the oldest reconstituted vial for the pool. Alternatively, prepare the necessary volumes of diluted conjugate and pool before using the assay.
2. Determine the number of wells to be used.
3. Prepare the enzyme conjugate solution by diluting the 180X enzyme conjugate stock 1 to 180 in the EIA Buffer provided. Mix the solution by inversion. Do not vortex. Example: for two eight well strips, add 25 µL of the 180X enzyme conjugate stock to 4475 µL of EIA Buffer.
4. Add 20 µL of sample, laboratory calibrators and CD controls to the appropriate wells in duplicate. Do not dilute CD's positive and negative controls.
5. Add 180 µL of diluted drug-enzyme conjugate to each well. Use 8-channel pipetter or 12-channel pipetter for rapid addition.
6. For manual runs, mix by gently shaking plate. A microplate shaker may be used.
7. Cover plate with plastic film or plate cover and incubate at room temperature for 45 minutes.
8. During the conjugate incubation, dilute concentrated wash buffer 10 fold with deionized water (i.e. 20 mL of concentrated wash buffer plus 180 mL of deionized water). Mix thoroughly. Diluted wash buffer is stable for 5 days at room temperature or 7 days at 2-8°C.
9. Once the incubation is complete, dump or aspirate the liquid from the wells. Tap the plate on a clean lint-free towel to remove any remaining liquid in the wells.
10. Wash each well with 300 µL of diluted wash buffer. Manual Wash: For manual wash procedures repeat for a total of 3 washings, invert and tap dry the plate following each step. After completing the last wash step wipe the bottom of the wells with a lint-free towel to remove any liquid on the outside of the wells. Automated

Wash: If an automated plate washer is used wash the plate for a total of 5 washings with 300 μ L of diluted wash buffer. It is important for the automated washer to conduct a final aspirate cycle to eliminate residual amounts of wash buffer. Residual amounts of buffer in the wells will affect assay performance. Note: DI water should never be used for the plate wash.

11. Add 150 μ L of the K-Blue Substrate to each well. For manual runs, use a multi-channel pipetter for best results.
12. Incubate at room temperature for 30 minutes. Gently shake immediately before measuring the absorbance.
13. Add 50 μ L of CD's Red Stop Solution to each well to stop enzyme reaction. Mix gently before measuring the absorbance. For automated systems a 10 second shake is sufficient. Measure the absorbance at a wavelength of 650 nm. Wells should be read within 2 hours of stopping the reaction. Note: When CD's Red Stop is used, it will result in a color ranging from a dark blue/purple to light pink based on the concentration of drug present. If 1 N HCl is preferred, use 50 μ L per well and read plate with a 450 nm filter. All QC data is generated without using a stopping reagent. Note: When acid stop is used, OD values will approximately double as compared to the OD values obtained with Red Stop.

Interpretation Of Results

Each laboratory should determine the cutoff level for their individual application. When possible, cutoff calibrators and/or standards should be prepared in the same matrix being tested.


Positive Result: Samples with an absorbance less than or equal to the laboratory's designated cutoff calibrator should be considered positive. All positive samples should be confirmed by a quantitative method such as GC/MS.

Negative Result: Samples with an absorbance greater than the laboratory's designated cutoff calibrator should be considered negative.

Qualitative QC Controls: The positive and negative controls provided in the kit are for QC purposes only. The sole purpose of these controls is to verify that the test kit is performing properly. The controls are not intended for use as cutoff calibrators. The positive control is spiked at a high concentration and its approximate level can be found on the label.

Note: The kit was designed for screening purposes only. It is recommended that all suspect samples be confirmed by a quantitative method such as GC/MS or HPLC.

Typical Standard Curve

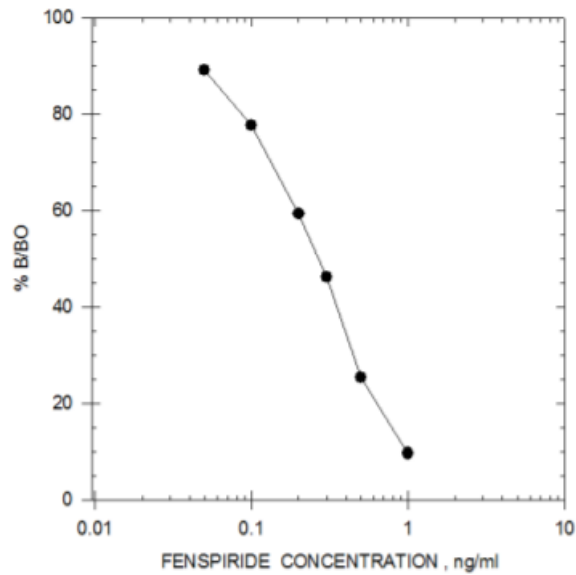

Sensitivity

Compound	I-50 in EIA Buffer
Fenspiride	0.09 ng/mL
Fentanyl	1.01 ng/mL
3-Methylfentanyl	2.3 ng/mL
Carfentanil	2.7 ng/mL
Sufentanil	21 ng/mL
Lofentanil	25 ng/mL
α -Methylfentanyl	286 ng/mL
Isoxsuprine	625 ng/mL
Promazine	1496 ng/mL

The term I-50 is used to define the sensitivity of the test. This number is derived from a standard curve generated with the drug in EIA Buffer. The drug concentration that shows 50% less color activity than the zero standard is considered to be the I-50.

Specificity

Compound	Compound Concentration (ng/mL)	Fenspiride Equivalents (ng/mL)	% Cross-Reactivity
Fenspiride	0.09	0.09	100%
Fentanyl	1.00	0.09	9%
3-Methylfentanyl	2.25	0.09	4%
Carfentanil	3.00	0.09	3%
Sufentanil	22.5	0.09	0.4%
Lofentanil	22.5	0.09	0.4%
α -Methylfentanyl	300	0.09	0.03%
Isoxsuprine	900	0.09	0.01%
Promazine	900	0.09	0.01%

Note: Fenspiride equivalents represents 50% B/B₀ assay displacement in EIA Buffer.

The compounds having cross-reactivity below 0.01% did not show any significant reaction up to 10 μ g/mL.

ALL THE FOLLOWING HAVE A CROSS-REACTIVITY <0.01%.

Acepromazine; Acetaminophen; Acetylsalicylic Acid; Amitriptyline; Amphetamine; Ascorbic Acid; Caffeine; Chlordiazepoxide; Chlorpromazine; Clenbuterol; Codeine; Cortisol; Cotinine; Dexamethasone; Dextromethorphan; Diclofenac; Dimethyl Sulfoxide; Dipyrone; Doxepin; E-amino-n-caproic Acid; Ephedrine; Erythromycin; Ethyl p-amino-benzoate; Fenopropfen; Flunixin; Folic Acid; Folinic Acid; Furosemide;

Gemfibrozil; Gentisic Acid; Glipizide; Glutethimide; Glycopyrrolate; Heparin; Hippuric Acid; Hordenine; Hydrocortisone; Hydromorphone; Ibuprofen; Imipramine; Isoxsuprine; L-Glutamic Acid; Lidocaine; Meperidine; Metaproterenol; Methadone; Methamphetamine; Methaqualone; Methocarbamol; Methylene Blue; Methylprednisolone; Nalorphine; Naproxen; Niacinamide; Nicotine; Nortriptyline; Orphenadrine; Oxycodone; Oxyphenbutazone; Penicillin G-Potassium; Penicillin G-Procaïne; Pentoxifylline; Phencyclidine (PCP); Phenothiazine; Phenylbutazone; Polyethylene Glycol; Prednisolone; Primadone; Procainamide; Procaine; Promazine; Pseudoephedrine; Pyrantel; Pyrilamine; Pyrimethamine; Quinidine; Quinine; Salbutamol; Salicylamide; Salicylic Acid; Theophylline; Thiamine; Trimethoprim, Trimpramine, Uric Acid.

Precautions

1. DO NOT use kits or components beyond expiration date.
2. DO NOT mix conjugates and plates from different kit lots.
3. DO NOT pipette reagents by mouth.
4. Pour K-Blue Substrate out of the bottle into a clean reservoir. To prevent contamination of the substrate, DO NOT pipette out of the bottle.
5. All specimens should be considered potentially infectious. Exercise proper handling precautions.
6. Keep plate covered except when adding reagents, washing or reading.
7. Kit components should be refrigerated at all times when not in use.
8. Keep the controls frozen if storing longer than 10 days. Avoid repeated freeze-thaw cycles. Note: Some kits require controls to be stored frozen immediately upon receipt. Reference kit label for details.
9. Use aseptic technique when opening and removing reagents from vials and bottles.
10. DO NOT smoke, eat or drink in areas where specimens or reagents are being handled.
11. Do not substitute DI water for the wash step of this protocol. Use only CD's wash buffer.
12. Do not use Sodium Azide with samples, standards and/or calibrators.
13. Do not reuse wells, they are for one use only.