



User's Manual

# Human Prostatic Acid Phosphatase (PAP) ELISA Kit



IVDIA1003-FA



96T





This product is for research use only and is not intended for diagnostic use.

For illustrative purposes only. To perform the assay the instructions for use provided with the kit have to be used.

---

Creative Diagnostics

 Address: 45-1 Ramsey Road, Shirley, NY 11967, USA

 Tel: 1-631-624-4882 (USA) 44-161-818-6441 (Europe)  Fax: 1-631-938-8221

 Email: [info@creative-diagnostics.com](mailto:info@creative-diagnostics.com)  Web: [www.creative-diagnostics.com](http://www.creative-diagnostics.com)

---

## PRODUCT INFORMATION

### Intended Use

Human Prostatic Acid Phosphatase (PAP) ELISA Kit is a quantitative solid phase enzyme-linked immunosorbent assay. This test provides quantitative measurement of human prostatic acid phosphatase (PAP) in serum.

### General Description

The phosphatases include two main types. The alkaline phosphatase has a pH optimum of about 9 while acid phosphatase (AcP) has its optimal activity at a pH of about 5. The AcP, first demonstrated in the urine in 1925, was found to be much more prevalent in male than in female urine. Prostatic acid phosphatase (PAP) enzyme activity was found to be localized in organs of the male genital tract. The presence of PAP in sera of patients with prostatic cancer was first demonstrated by Gutman. It was shown that PAP concentration was elevated in many men with primary prostatic carcinoma and metastatic lesions of the prostate.

The release of PAP into the blood stream by primary or metastatic carcinoma of the prostate may provide the clinician with a means of following remission or relapse of a prostatic malignancy. Therefore, sensitive and accurate measurement of serum PAP is essential in monitoring the effectiveness of various therapeutic treatments.

Because moderate elevation of serum PAP is observed to accompany many nonprostatic diseases, a specific assay for determination of PAP is established. The current methodologies for determination of PAP activity using the colorimetric methods and immunochemical methods have many limitations.

### Principles of Testing

Human Prostatic Acid Phosphatase (PAP) ELISA Kit is a solid phase enzyme linked immunosorbent assay (ELISA). The wells are coated with specific anti-PAP antibodies. The Samples, Standards and Controls are incubated in the wells and bound with enzyme conjugate which is a mixture of anti-PAP antibodies with different affinity toward epitopes of PAP molecule and chemically conjugated with horseradish peroxidase. Unbound enzyme conjugate is washed off. The amount of bound peroxidase is proportional to the concentration of the PAP present in the Samples, Standards and Controls. Upon addition of the TMB substrate, the intensity of color developed is proportional to the concentration of PAP.

### Reagents And Materials Provided

1. Micro-well Strips: Rabbit anti-PAP antibodies coated wells, 96 wells.
2. Enzyme Conjugate (11 ml): Mouse IgG(anti-PAP) conjugated to horseradish peroxidase
3. Sample Diluent or Zero Standard (11 mL).
4. Reference Standard Set (0.75 mL/each): Calibrated to 1, 3, 5, 15, and 30 ng/mL in BSA-containing diluent.
5. TMB Solution (11 mL): Buffer solution containing hydrogen peroxide and TMB
6. Washing buffer concentrate 100x (10mL): Prepare working washing solution by adding 10 mL washing buffer concentrate into 990 mL distilled water.

7. Stop Solution: 2 N HCl.
8. Well holder: for securing wells.
9. Package Insert.

## Materials Required But Not Supplied

1. Micro-well reader.
2. Disposal tips and pipettor of 25  $\mu$ L and 100  $\mu$ L.
3. Quality controls such as Bio-Rad Lyphocheck were recommended

## Storage

1. Store the kits at 2-8°C in a refrigerator.
2. Keep micro-wells in a dry bag with desiccants.
3. The reagents are stable until expiration of the kit. TMB Solution should be colorless; if the solution turns blue, it must be replaced.

## Specimen Collection And Preparation

Collect blood by venipuncture, allow to clot, and separate the serum by centrifugation at room temperature. Serum is required for PAP assay, and do not add sodium azide as preservative. PAP is unstable at 4°C and 20°C. Samples should not be stored at room temperature or 4°C for more than 24 hours. Serum samples are recommended to be frozen for longer storage. Avoid repeated freezing and thawing of serum samples.

## Assay Procedure

### Preparation for Assay

1. Bring all reagents and samples to room temperature (24 $\pm$ 3°C) and shake gently before beginning the test. Have all reagents and samples ready before the start of the assay. Once the test is begun it must be performed without any interruption to get the most reliable and consistent results.
2. Use new disposable tips for each specimen.

### Assay Procedure

1. Secure the desired number of coated wells in the holder.
2. Dispense 25  $\mu$ L of references, controls or serum samples into the appropriate wells.
3. Dispense 100  $\mu$ L of enzyme conjugate into wells. Incubate for 30 minutes at room temperature.
4. Remove incubation mixture and rinse the wells 5 times with washing buffer (300  $\mu$ L/well/each rinse).
5. Dispense 100  $\mu$ L of TMB Solution into each well and incubate for 15 minutes at room temperature.
6. Stop reaction by adding 50  $\mu$ L of stop solution into each well.
7. Read O.D. at 450 nm with a microwell reader in 5 minutes.

### Note:

1. Wash the microwells and remove washing buffer thoroughly.
2. Pipet all reagents and samples into bottom of the well. Vortex-mixing or shaking is not required.
3. Absorbance is a function of the time and temperature of incubations. It is recommended to have reagents, samples and needed wells ready for ensure the equal elapsed time for each pipetting without interruption.
4. For the same reason run no more than 20 patient samples with a set of reference standards in duplicate for each assay.
5. If a serum specimen contains greater than 30 ng/ml of PAP, the sample must be diluted with sample diluent and re-assayed as described in the assay procedure.


### Quality Control

Good laboratory practices include the use of control specimens to ensure that all reagents and protocols are performing properly. Assay do not include serum controls. Bio-Rad Lyphochek series control serum is recommended.

### Calculation

1. Plot the concentration (X) of each reference standards against its absorbance (Y) on a full logarithmic graph paper.
2. Obtain the PAP value of patient by reference to the standard curve as follows: (These data are for demonstration purpose only and must not be used in place of data generated for each assay).

| Well No. | Description (ng/mL) | Absorbance 450 nm | PAP (ng/mL) |
|----------|---------------------|-------------------|-------------|
| A1       | 0                   | 0.000             |             |
| B1       | (Blank)             | 0.000             |             |
| A2       | 1                   | 0.102             |             |
| B2       |                     | 0.089             |             |
| A3       | 3                   | 0.215             |             |
| B3       |                     | 0.241             |             |
| A4       | 5                   | 0.371             |             |
| B4       |                     | 0.383             |             |
| A5       | 15                  | 1.019             |             |
| B5       |                     | 1.039             |             |
| A6       | 30                  | 2.054             |             |
| B6       |                     | 2.060             |             |
| A7       | PATIENT A           | 1.653             | 23.19       |
| B7       |                     | 1.455             |             |

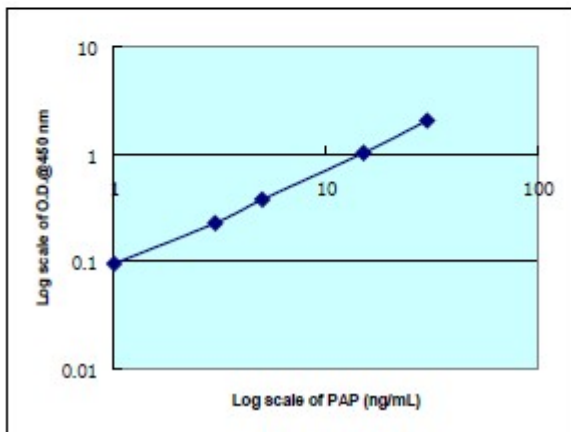

## Reference Values

1. Normal serum concentrations of prostatic acid phosphatase in men have been reported to range from zero to 5 ng/mL. 2. Elevated serum concentrations have been observed frequently in patients with prostatic carcinoma and metastatic lesions of the prosta

## Specificity

The test is specific for human PAP. Cross-reactivity with other serum proteins, human hormones and the tumor makers were not found.

## Recovery

Recovery studies were performed by mixing an aliquot of pooled serum and PAP Standards. The PAP values were measured and percentage of recovery determined.

## Precautions

1. It is designed for in Vitro Diagnostic Use only.
2. The Components in the kits are intended for usage as an integral unit. The components from different lots should not be mixed, and not be used beyond expiration date.
3. The material should be used in a designated work area, the bench surface should be cleaned with detergent and the contaminated materials should be disposed properly.
4. Some components have been tested using FDA-approved methods and has been found negative for antibody to human immuno-deficiency virus (HIV-I, HIV-II), antibody to Hepatitis C and Hepatitis B surface antigen (HBsAg). No known test method can offer total assurance that HIV-I, HIV-II, Hepatitis B & C virus or other infectious agents are absent. Handle these reagents as if they were potentially infectious.
5. Avoid microbial contamination of reagents when removing aliquots from the vials.

## Limitations

1. For diagnostic purpose, the PAP values should be used as an adjunct to other data available to the physician.
2. The kit is designed to avoid "hook effect" up to 1000 ng/mL.
3. Samples with PAP level above 30 ng/mL should be diluted to obtain accurate v

## References

1. Gorman, L. and Statland, B.E. Clinical usefulness of alkaline phosphatase Isoenzyme determinations. Clin. Biochem., 10:171, 1977.
2. Kutscher, W. and Wolbergs H.. Prostatic Phosphatase I. Z. Physiol Chem. 236:237, 1935.
3. Kutscher, W. and Warner, A. P