



**User's Manual** 

# **BSA ELISA Kit**

REF

**DEIA-PX0002** 



96T



This product is for research use only and is not intended for diagnostic use.

For illustrative purposes only. To perform the assay the instructions for use provided with the kit have to be used.

#### **Creative Diagnostics**

Address: 45-1 Ramsey Road, Shirley, NY 11967, USA

Tel: 1-631-624-4882 (USA) 44-161-818-6441 (Europe)

Fax: 1-631-938-8221

Email: info@creative-diagnostics.com

Web: www.creative-diagnostics.com

## PRODUCT INFORMATION

#### **Intended Use**

The BSA ELISA Detection Kit is designed for the quantitative detection of residual BSA content in intermediates, semi-finished products and finished products of various biological products

# **Principles of Testing**

The kit can quantitative detect residual BSA content in test samples by using a double-antibody sandwich method. Coat the 96-well plate with capture antibody to prepare immobilized antibodies. Add the standard and test samples, and then horseradish peroxidase (HRP) conjugated antibody to form an immobilized antibody - BSA - enzyme-labeled antibody sandwich conjugate. Wash the plate after reaction, and add the substrate for color development reaction. The substrate will turn blue under HRP catalysis, and will finally turn yellow under the action of the stop solution. Determine the optical density (OD) value at 450 nm, and calculate the BSA content in the test sample using the standard curve.

# **Reagents And Materials Provided**

- BSA Coated Plate, 8 wells x 12 strips, Ready-to-use 1.
- 2. BSA standard (standard), 100 µL x 1 vial (5 µg/mL), Operate as per the recommended dilution procedure
- 3. 1250x Anti-BSA (1250x enzyme-labeled antibody), 30 μL x 1 vial, Operate as per the recommended dilution procedure
- 20x Buffer, 30 mLx 1 bottle, Operate as per the recommended dilution procedure 4.
- 5. Color Reagent A, 8 mLx 1 vial, Ready-to-use
- 6. Color Reagent B, 8 mLx 1 vial, Ready-to-use
- 7. Stop Solution, 15 mLx 1 bottle, Ready-to-use
- 8. Plate Sealer, 3 pieces, Ready-to-use
- 9. Instructions for Use, 1 copy

Notes: Reagents should be stored at 2 ~ 8°C; Color Reagents A and B should be protected from light during storage

### Materials Required But Not Supplied

- Plate reader, thermostate plate shaker or thermostate incubator, plate washer (should not be cross-used in other rojects).
- 2. High-precision pipette and disposable tips  $(0.5 \sim 10 \,\mu\text{L}, 10 \sim 100 \,\mu\text{L}, 30 \sim 300 \,\mu\text{L}, \text{ and } 100 \sim 1000 \,\mu\text{L})$ .
- 3. Deionized water, blotting paper, EP tubes

#### **Storage**

The shelf life of unopened kits is 12 months when stored at  $2 \sim 8^{\circ}$ C.

Tel: 1-631-624-4882 (USA)

Tel: 44-161-818-6441 (Europe)




Fax: 1-631-938-8221

Email: info@creative-diagnostics.com

# **Reagent Preparation**

- Avoid placing electronic balances or weighing BSA powder in the laboratory as much as possible. Clean the benchtop before experiments.
- 2. Preparation of 1x Buffer: Take 1 portion of Buffer (20x), and add 19 times the volume of deionized water to prepare the working concentration buffer (1x). If there are crystals in the Buffer (20x), shake gently at room temperature or in a 37°C water bath, and dilute after the crystals are completely dissolved. Unused Buffer (20x) should be stored at 2 ~ 8°C. **Notes:** Use clean, dedicated containers to prepare the Buffer as needed, and do not share solutions, containers, or consumables with other projects.
- Preparation of standard: Dilute the standard to 50 ng/mL with 1x Buffer, and prepare the standard by 2-fold dilution, as shown in the figure above.

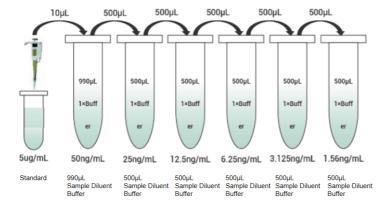

- Preparation of 1x enzyme-labeled antibody: Take out the 1250x enzyme-labeled antibody from the freezer and place it on an ice box. Dilute the 1250x enzyme-labeled antibody with 1x Buffer by 1250 folds based on the volume required (100 μL/well) to prepare 1x enzyme-labeled antibody. Do not leave the 1250x enzymelabeled antibody at room temperature for a long time. It is preferred to take the antibody out when needed and perform operations on an icebox. The 1x enzyme-labeled antibody should be prepared freshly before use.
- Preparation of substrate solution: Mix Color Reagents A and B at equal volume at 10 minutes before use, and the operation should be performed at dark environment. Make sure that the substrate solution is not contaminated. Do not use if the substrate solution turns blue after mixing.

# **Assay Procedure**

- Equilibrate the temperature of each component in the kit to room temperature for 30 minutes. Take out required strip plates from aluminum foil bags already equilibrated to room temperature, and label the strip plate sequence with a marker. Seal remaining strip plates with a plate sealer, put them back to the aluminum foil bag, then seal the bag, and store at  $2 \sim 8^{\circ}$ C.
- Sample incubation: Set the standard wells, blank wells, and test sample wells, respectively. Add standards 2. at different concentrations (in sequence), 1x Buffer, and test sample to standard wells, blank wells, and test sample wells, respectively (100 µL/well; replicate well is recommended to avoid drifts), seal the plate with a plate sealer, and incubate at 37°C for 1 hour.
- Plate washing: Discard liquid in the wells; wash the plate for 3 times with 1x Buffer (300 μL/well), and pat dry the residual liquid in wells on a blotting paper. (If the plate is washed manually, avoid the pipette tip touching the inner wall of each well when adding 1x Buffer for plate washing. After adding Buffer each time,

Tel: 1-631-624-4882 (USA)

Tel: 44-161-818-6441 (Europe)



Fax: 1-631-938-8221


allow to stand for 1 minute and shake gently. When patting for drying, use a new blotting paper each time, or dry the plate in a clean area of the paper)

- Incubation of enzyme-labeled antibody: Add 1x enzyme-labeled antibody at 100 μL/well, seal the plate with a plate sealer, and incubate at 37°C for 1 hour.
- 5. Plate washing: Same as Step 3; wash the plate for 5 times.
- Color development: Add the pre-prepared substrate solution into the plate at 100 µL/well, seal the plate with a plate sealer, and allow the solution to stand and incubate at 37°C for 20 minutes while being protected from light.
- 7. Reaction termination: Add stop solution at 100 μL/well.
- Reading: Measure the OD values at 450 nm and 630 nm with a plate reader. The measurement should be completed within 20 minutes after reaction termination.

#### Calculation

Plot the standard curve with the OD value (OD450 nm - OD630 nm) of the BSA standard as the dependent variable Y and the standard concentration (ng/mL) as the independent variable X. The 4-parameter logistic fitting equation is recommended:  $Y = ((A - D)/(1 + (x/C)^B)) + D$ . Substitute the OD value of test sample (OD450 nm - OD630 nm) to the equation of the standard curve to calculate the BSA content in test sample.

# **Typical Standard Curve**

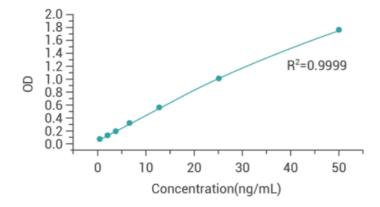

Standard concentration (ng/mL)	OD value (1)	OD value (2)	Mean value
50	1.775	1.803	1.789
25	1.091	0.995	1.043
12.5	0.571	0.591	0.581
6.25	0.324	0.308	0.316
3.125	0.193	0.179	0.186
1.56	0.132	0.127	0.130
0	0.072	0.088	0.080

#### **Precautions**

The containers used for BSA assay should be dedicated, so as to prevent BSA contamination in the laboratory.

Tel: 1-631-624-4882 (USA)

Tel: 44-161-818-6441 (Europe)



Fax: 1-631-938-8221

Email: info@creative-diagnostics.com

- 2. This kit is for in vitro detection only, and may not be used for clinical diagnosis.
- 3. The kit must be used within the shelf life. Mix well each component before use.
- 4. To ensure accurate assay, please strictly follow the instructions for use, use only the reagents provided in the kit, and do not use self-prepared reagents or reagents from other kits.
- 5. Please timely replace reagent troughs and pipette tips when loading different samples and performing different steps, so as to avoid cross contamination.
- The final assay results are closely related to reagent effectiveness, the operations of analysts, and the test 6. environment.
- 7. Our company is only responsible for the kits themselves, and will not be responsible for the sample consumption caused by kits during use. Users should fully consider the possible sample consumption before operation, and should reserve sufficient sample size.
- The stop solution in this kit is acidic, thus extra attention should be paid during operation.
- All biological samples may cause potential biosafety risks, therefore, users must strictly follow local laws and regulations when handling and disposing of the samples.

Tel: 1-631-624-4882 (USA)

Tel: 44-161-818-6441 (Europe)



**Fax:** 1-631-938-8221

Email: info@creative-diagnostics.com