



## User's Manual

# Thiostatin ELISA Kit



DEIA3846



96T



This product is for research use only and is not intended for diagnostic use.

For illustrative purposes only. To perform the assay the instructions for use provided with the kit have to be used.

---

### Creative Diagnostics

 Address: 45-1 Ramsey Road, Shirley, NY 11967, USA

 Tel: 1-631-624-4882 (USA) 44-161-818-6441 (Europe)  Fax: 1-631-938-8221

 Email: [info@creative-diagnostics.com](mailto:info@creative-diagnostics.com)  Web: [www.creative-diagnostics.com](http://www.creative-diagnostics.com)

---

## PRODUCT INFORMATION

### Intended Use


The Thiostatin ELISA Kit is a highly sensitive two-site enzyme linked immunoassay (ELISA) for measuring THIOSTATIN in biological samples of rats.

### General Description

Thiostatin is a major protease inhibitor in serum and an acute phase protein that increases significantly in concentration in the rat as a result of inflammation. The major pathophysiological role for rat Thiostatin has yet to be conclusively defined.

### Principles of Testing

The principle of the double antibody sandwich ELISA is represented in following figure. In this assay the THIOSTATIN present in samples reacts with the anti-THIOSTATIN antibodies which have been adsorbed to the surface of polystyrene microtitre wells. After the removal of unbound proteins by washing, anti-Thiostatin antibodies conjugated with horseradish peroxidase (HRP), are added. These enzyme-labeled antibodies form complexes with the previously bound Thiostatin. Following another washing step, the enzyme bound to the immunosorbent is assayed by the addition of a chromogenic substrate, 3,3',5,5'-tetramethylbenzidine (TMB). The quantity of bound enzyme varies directly with the concentration of Thiostatin in the sample tested; thus, the absorbance, at 450 nm, is a measure of the concentration of Thiostatin in the test sample. The quantity of Thiostatin in the test sample can be interpolated from the standard curve constructed from the standards, and corrected for sample dilution.

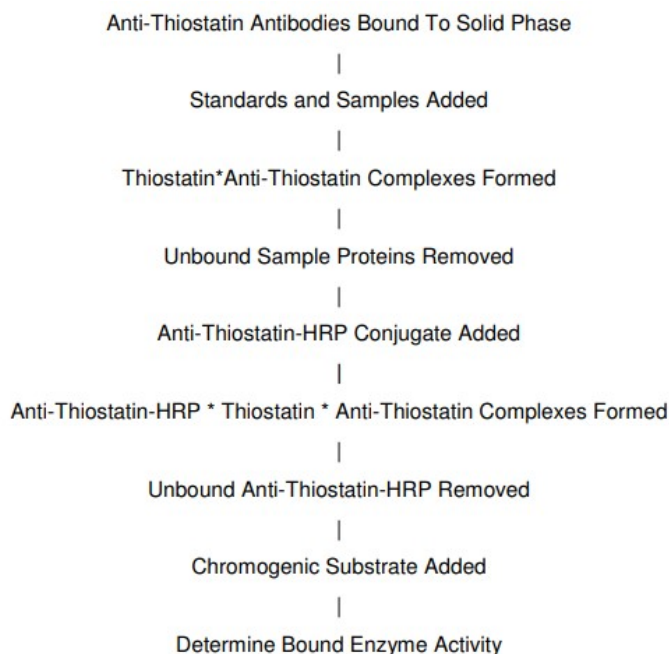

### Reagents And Materials Provided

1. Diluent Concentrate (Running Buffer): One bottle containing 50 ml of a 5× concentrated diluent running buffer.
2. Wash Solution Concentrate: One bottle containing 50 ml of a 20× concentrated wash solution.
3. Enzyme-Antibody Conjugate 100×: One vial containing 150 µl of affinity purified anti-Rat THIOSTATIN antibody conjugated with horseradish peroxidase in a stabilizing buffer.
4. Chromogen-Substrate Solution: One vial containing 12 mL of 3,3',5,5'-tetramethybenzidine (TMB) and hydrogen peroxide in citric acid buffer at pH 3.3.
5. Stop Solution: One vial containing 12 ml 0.3 M sulfuric acid. Warning: Avoid contact with skin.
6. Anti-Rat Thiostatin Elisa Micro Plate: Twelve removable eight (8) well micro well strips in well holder frame. Each well is coated with affinity purified anti-Rat Thiostatin.
7. Rat Thiostatin Calibrator: One vial containing Rat Thiostatin calibrator. For In Vitro Use Only

## Materials Required But Not Supplied

1. Precision pipette (2 µl to 200 µl) for making and dispensing dilutions
2. Test tubes
3. Microtitre washer/aspirator
4. Distilled or Deionized H<sub>2</sub>O
5. Microtitre Plate reader
6. Assorted glassware for the preparation of reagents and buffer solutions
7. Timer
8. Vortex mixer

## Storage

The expiration date for the package is stated on the box label.

1. Diluent: The 5× Diluent Concentrate is stable until the expiration date. The 1× working solution is stable for at least one week from the date of preparation. Both solutions should be stored at 4-8°C.
2. Wash Solution: The 20× Wash Solution Concentrate is stable until the expiration date. The 1× working solution is stable for at least one week from the date of preparation. Both solutions can be stored at room temperature (16-25°C) or at 4-8°C.
3. Enzyme-Antibody Conjugate: Undiluted horseradish peroxidase anti-Thiostatin conjugate should be stored at 4-8°C and diluted immediately prior to use. The working conjugate solution is stable for up to 8 hours.
4. Chromogen-Substrate Solution: The Substrate Solution should be stored at 4-8°C and is stable until the expiration date.
5. Stop Solution: The Stop Solution should be stored at 4-8°C and is stable until the expiration date.
6. Anti- Rat Thiostatin Elisa Micro Plate: Anti-Rat Thiostatin coated wells are stable until the expiration date, and should be stored at 4-8°C in sealed foil pouch with desiccant pack.
7. Rat Thiostatin Calibrator: Long Term Storage: Upon receipt, aliquot the calibrator and store them frozen. They will be stable until expiration date. Short Term Storage: the calibrator is stable for up to 14 days at 4°C. The working standard solutions should be prepared immediately prior to use and are stable for up to 8

hours.

## Specimen Collection And Preparation

Blood should be collected by venipuncture. The serum should be separated from the cells after clot formation by centrifugation. For plasma samples, blood should be collected into a container with an anticoagulant and then centrifuged. Care should be taken to minimize hemolysis, excessive hemolysis can impact your results. Assay immediately or aliquot and store samples at -20°C. Avoid repeated freeze-thaw cycles.

1. **Precautions:** For any sample that might contain pathogens, care must be taken to prevent contact with open wounds.
2. **Additives and Preservatives:** No additives or preservatives are necessary to maintain the integrity of the specimen. Avoid azide contamination.
3. **Known interfering substances** Azide and thimerosal at concentrations higher than 0.1% inhibits the enzyme reaction.

## Dilution of Samples

The assay for quantification of Thiostatin in samples requires that each test sample be diluted before use. For a single step determination a dilution of 1/20,000 is appropriate for most serum/plasma samples. For absolute quantification, samples that yield results outside the range of the standard curve, a lesser or greater dilution might be required. If unsure of sample level, a serial dilution with one or two representative samples before running the entire plate is highly recommended.

To prepare a 1/20,000 dilution of sample, transfer 5 µl of sample to 495 µl of 1× diluent. This gives you a 1/100 dilution. Next, add 5 µl of your 1/100 diluted sample to 995 µl of 1× diluent. You now have a 1/20,000 dilution of your sample. Mix thoroughly at each stage.

## Reagent Preparation

1. **Diluent Concentrate:** The Diluent Solution supplied is a 5× Concentrate and must be diluted 1/5 with distilled or deionized water (1 part buffer concentrate, 4 parts dH<sub>2</sub>O).
2. **Wash Solution Concentrate:** The Wash Solution supplied is a 20× Concentrate and must be diluted 1/20 with distilled or deionized water (1 part buffer concentrate, 19 parts dH<sub>2</sub>O). Crystal formation in the concentrate is not uncommon when storage temperatures are low. Warming of the concentrate to 30-35°C before dilution can dissolve crystals.
3. **Enzyme-Antibody Conjugate:** Calculate the required amount of working conjugate solution for each microtitre plate test strip by adding 10 µl Enzyme-Antibody Conjugate to 990 µl of 1× Diluent for each test strip to be used for testing. Mix uniformly, but gently. Avoid foaming.
4. **Chromogen-Substrate Solution:** Ready to use as supplied.
5. **Stop Solution:** Ready to use as supplied.
6. **Anti-Rat Thiostatin Elisa Micro Plate:** Ready to use as supplied. Unseal Microtiter Pouch and remove plate from pouch. Remove all strips and wells that will not be used in the assay and place back in pouch and re-seal along with desiccant.
7. **Rat Thiostatin Calibrator:** The calibrator is now at a concentration of 800 ng/ml. Rat Thiostatin standards need to be prepared immediately prior to use (see chart below). Mix well between each step. Avoid foaming.

Standard	ng/mL	Volume added to 1X Diluent	Volume of 1X Diluent
7	400	300 µl Thiostatin Calibrator	300 µl
6	200	300 µl standard 7	300 µl
5	100	300 µl standard 6	300 µl
4	50	300 µl standard 5	300 µl
3	25	300 µl standard 4	300 µl
2	12.5	300 µl standard 3	300 µl
1	6.25	300 µl standard 2	300 µl
0	0		500 µl

## Assay Procedure

1. Bring all reagents to room temperature before use.
  2. Pipette 100 µl of  
Standard 0 (0.0 ng/ml) in duplicate  
Standard 1 (6.25 ng/ml) in duplicate  
Standard 2 (12.5 ng/ml) in duplicate  
Standard 3 (25 ng/ml) in duplicate  
Standard 4 (50 ng/ml) in duplicate  
Standard 5 (100 ng/ml) in duplicate  
Standard 6 (200 ng/ml) in duplicate  
Standard 7 (400 ng/ml) in duplicate
  3. Pipette 100 µl of sample (in duplicate) into pre designated wells.
  4. Incubate the micro titer plate at room temperature for sixty ( $60 \pm 2$ ) minutes. Keep plate covered and level during incubation.
  5. Following incubation, aspirate the contents of the wells.
  6. Completely fill each well with appropriately diluted Wash Solution and aspirate. Repeat three times, for a total of four washes. If washing manually: completely fill wells with wash buffer, invert the plate then pour/shake out the contents in a waste container. Follow this by sharply striking the wells on absorbent paper to remove residual buffer. Repeat 3 times for a total of four washes.
  7. Pipette 100 µl of appropriately diluted Enzyme-Antibody Conjugate to each well. Incubate at room temperature for sixty ( $60 \pm 2$ ) minutes. Keep plate covered in the dark and level during incubation.
  8. Wash and blot the wells as described in Steps 5/6.
  9. Pipette 100 µl of TMB Substrate Solution into each well.
  10. Incubate in the dark at room temperature for precisely ten (10) minutes.
  11. After ten minutes, add 100 µl of Stop Solution to each well.
  12. Determine the absorbance (450 nm) of the contents of each well. Calibrate the plate reader to air. H.
- The absorbance of the final reaction mixture can be measured up to 2 hours after the addition of the Stop Solution. However, good laboratory practice dictates that the measurement be made as soon as possible.

## Calculation

1. Subtract the average background value from the test values for each sample.
2. Using the results observed for the standards construct a Standard Curve. The appropriate curve fit is that of a four-parameter logistics curve. A second order polynomial (quadratic) or other curve fits may also be used.
3. Interpolate test sample values from standard curve. Correct for sera dilution factor to arrive at the Thiostatin concentration in original samples.

## Limitations

1. Reliable and reproducible results will be obtained when the assay procedure is carried out with a complete understanding of the information contained in the package insert instructions and with adherence to good laboratory practice.
2. Factors that might affect the performance of the assay include proper instrument function, cleanliness of glassware, quality of distilled or deionized water, and accuracy of reagent and sample pipettings, washing technique, incubation time or temperature.
3. Do not mix or substitute reagents with those from other lots or sources.

